

ADVANCED

SPIROBANK II ADVANCEDTM

Handheld, Stand-alone and
PC-based Spirometer, with
Oximetry Option

A complete Spirometer and Oximeter,
designed for all respiratory therapists.



MAIN features



REAL-TIME TEST

Spirometry: FVC, VC, IVC, MVV, PRE/POST
Bronchodilator comparison
Oximetry (optional): Spot test (SpO2, BPM)



CARRY EVERYWHERE

High resolution backlight display, long battery life, large internal storage, carrying case included



COMPLIANCE ATS/ERS 2019

And other Standards including ISO 26782 (for Spirometry), ISO 23747 (for PEF), ISO 80601-2-61 (for Oximetry), and more. CE0476, FDA 510 (k)



PC CONNECTION AVAILABLE

Real-time test on PC screen, connect with your EHR/EMR, back-up internal memory and more, via USB and Bluetooth



DISTINCTIVE features



PREDICTED SETS & VALUES

Large Selection, including comparison %Pred, Z-score and LLN. Include GLI in PC-mode



OXIMETRY OPTION

Spot test in Real-Time



EHR/EMR CONNECTIVITY

Via PC, integration with patient database on your EHR/EMR (in HL7, GDT)



COVID-19 PREVENTION

Complete Disposable Set with Antiviral filter. Bluetooth connection to test at safety distance

Always INCLUDED


- ✦ Carrying case
- ✦ USB cable
- ✦ Noseclip
- ✦ PC Software license

With Oximetry Option:

- ✦ Finger Probe

Compatible SOFTWARE

winspiroPRO

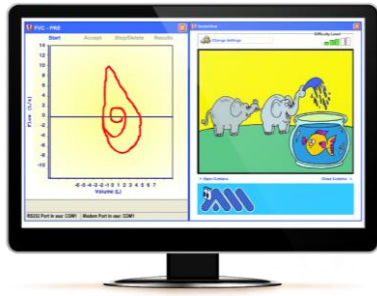

Pediatric Incentive (PATENTED) to improve patient compliance during the test.

Acceptability Messages, Test interpretation and Quality Control Grade according to the latest **Spirometry Standards**

MAIN FEATURES

Windows-based solution for Spirometry, Oximetry and Telemedicine.


Wide range of predicted sets and values, including **GLI Predicted sets, LLN and Z-score.**

Embedded **EHR/EMR connectivity.**


NET VERSION available, share one database between different PC workstations.

MEDICAL REPORT

Specialized and **customizable printout**

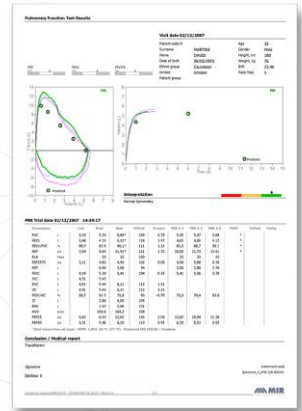

spiro Connect

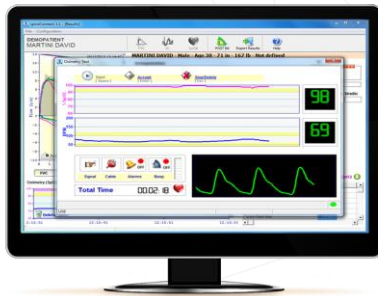

MAIN FEATURES

Windows-based solution, **direct integration** with your EHR/EMR.

Real time test include **Spirometry and Oximetry**
Standardized communication in **HL7 or Exchange Protocol.**

Select patient info directly from your own **EHR/EMR**

Spirometry test: FVC-Pre, FVC-Post, VC-Pre
Oximetry Test: SpO2 (%), Pulse (BPM)

GO-TO-MARKET TOOLKIT

Software Development Kit available for System Integrators and App Developers.
OEM service available for Spirometry and Oximetry.



Learn more about available SDK and OEM



Compatible TURBINES

flowMIR™
Disposable Turbine



Reusable Turbine



Mouthpiece

Included Disposable

Required, Not Included

Turbine Disinfection

Not required

Required

Turbine Calibration

Not required

Required

Packaging

Individually sealed: 60 or 10 units / box

1 unit in Carton box

Antiviral Filter

Available Disposable

Required Disposable

VISIT WEBSITE



SCIENTIFIC PUBLICATIONS



Also available in **MORE CONFIGURATIONS**



Technical Specification

Spirobank II Advanced

Spirobank II Smart

Spirobank II Basic

	Spirobank II Advanced	Spirobank II Smart	Spirobank II Basic
TYPE OF SPIROMETER	StandAlone + PC, with Oximetry Option	StandAlone + PC + App, with Oximetry Option	StandAlone + PC
COMPATIBLE TURBINES	flowMIR™ Disposable Turbine, Reusable Turbine Flowmeter	flowMIR™ Disposable Turbine, Reusable Turbine Flowmeter	flowMIR™ Disposable Turbine, Reusable Turbine Flowmeter
COMPATIBLE SOFTWARES	Winspiro PRO, spiro Connect	MIR Spiro App, Winspiro PRO, spiro Connect	Winspiro PRO, spiro Connect
EXTERNAL CONTROL	<p>Real-Time test on PC screen, connect with your EHR/EMR, back-up internal memory, and much more</p> <p>Connect to your PC via USB and Bluetooth 2.0</p>	<p>Real-Time test on Tablet screen and PC screen, connect with your EHR/EMR, back-up internal memory, and much more</p> <p>Connect to your PC via USB (no Bluetooth)</p> <p>Connect to your Tablet via Bluetooth Smart BLE 4.0</p>	<p>Real-Time test on PC screen, connect with your EHR/EMR, back-up internal memory and much more</p> <p>Connect to your PC via USB</p>
EHR CONNECTIVITY	Via PC, integration with patient database on your EHR/EMR (in HL7, GDT)	Via PC: integration with patient database on your EHR/EMR (in HL7, GDT) Via APP: transfer data to a remote server in HL7 standards	Via PC, integration with patient database on your EHR/EMR (in HL7, GDT)
MEASURED PARAMETERS	<p>Spirometry: FVC, VC, IVC, MVV, PRE-POST Bronchodilator comparison Oximetry (optional): Spot test (SpO2, BPM)</p> <p>Spirometry: FVC, FEV1, FEV1/FVC%, DTPEF, FEV 0.5, FEV0.5/FVC%, FEV0.75, FEV0.75/FVC%, FEV2, FEV2/FVC%, FEV3, FEV3/FVC%, FEV6, FEV1/FEV6%, PEF, FEF25%, FEF50%, FEF75%, FEF25-75%, FEF75-85%, FET, Vext, ELA, EVOL, FIVC, FIV1, PIF, FIV1/FIVC%, FIF25%, FIF50%, FIF75%, R50, PIF, IRV, VC, IVC, IC, ERV, FEV1/VC%, TV, VE, RR, ti, te, ti/t-tot, TV/ti, MVV</p> <p>Oximetry (Optional): SpO2% (min, max, average), BPM (min, max, average), Test duration, % Bradycardia Duration (<40 BPM), % Tachycardia Duration (>120 BPM), % of Time with SpO2 ≤ 90% (T90%, T89%)</p>	<p>Spirometry: FVC, VC, IVC, MVV, PRE-POST Bronchodilator comparison Oximetry (optional): Spot test (SpO2, BPM)</p> <p>Spirometry: FVC, FEV1, FEV1/FVC%, DTPEF, FEV 0.5, FEV0.5/FVC%, FEV0.75, FEV0.75/FVC%, FEV2, FEV2/FVC%, FEV3, FEV3/FVC%, FEV6, FEV1/FEV6%, PEF, FEF25, FEF50, FEF75, FEF25-75, FEF75-85%, FET, Vext, ELA, EVOL, FIVC, FIV1, PIF, FIV1/FIVC%, FIF25, FIF50, FIF75, R50, PIF, IRV, VC, IVC, IC, ERV, FEV1/VC%, TV, VE, RR, ti, tE, ti/t-tot, TV/ti, MVV</p> <p>Oximetry (Optional): SpO2% (min, max, average), BPM (min, max, average), Test duration, % Bradycardia Duration (<40 BPM), % Tachycardia Duration (>120 BPM), % of Time with SpO2 ≤ 90% (T90%, T89%)</p> <p>on MIR Spiro App: Spirometry: FVC, VC, PRE/POST Bronchodilator comparison Parameters: FVC, FEV1, FEV1%, PEF, FEF25-75, FET, Lung Age, VC, IVC. Oximetry (Optional): %SpO2 [Baseline, Min, Max, Mean], Pulse Rate [Baseline, Min, Max, Mean] Events.</p>	<p>Spirometry: FVC, VC, IVC, PRE/POST Bronchodilator comparison</p> <p>Spirometry: FVC, VC, IVC, IC, ERV, FEV1, FEV1%, PEF, FEF 25-75, FET, EVOL, ELA</p>

COMPARE ON WEBSITE



TECHNICAL datasheet

PRODUCT CODES - Spirobank II Advanced Configurations

911020E0 - Spirometer • 911020E1 - Spirometer with reusable turbine

911025E0 - Spirometer + Oximeter • 911025E1 - Spirometer + Oximeter with reusable turbine

Technical specification

Width	55 mm
Length	160 mm
Thickness	25 mm
Weight	140 g (battery pack included)

Turbine



Reusable turbine (code 910002)



Disposable turbine (code 910004)

Power supply

Rechargeable Lithium-Ion 3.7V,
1100 mAh

Current capacity

1100 mAh
~20-30 mA (during test)

Backup battery voltage

none

Batteries charger

voltage=5 V DC,
current=minimum 500 mA,
input current= 100VAC - 240 VAC
Connector: micro USB type B
compliant with EN 60601-1
50 hours

Autonomy

Connectivity

USB 2.0, Bluetooth® 2.1

Display

LCD monochrome, 160 × 80 pixel

Keyboard

membrane keyboard with 6 keys

Mouthpieces

Ø 30 mm (1.18 inch)

Type of electrical protection

Internally powered

Safety level for shock hazard

Type BF Apparatus

Conditions of use

Apparatus for continuous use

Conditions of storage

Temperature: MIN -20 °C,
MAX + 60 °C
Humidity: MIN 10% RH;
MAX 95%RH

Operating Conditions

Temperature: MIN + 10 °C,
MAX + 40 °C
Humidity: MIN 10% RH,
MAX 95%RH

Applied norms

Electrical Safety Standard
EN 60601-1
Electro Magnetic Compatibility
EN 60601-1-2

Spirometry

Flow sensor	bi-directional digital turbine
Flow range	±16L/s
Volume accuracy	±2.5% or 50 mL
Flow accuracy	±5% or 200 mL/s
Dynamic resistance	<0.5 cm H ₂ O/L/s
Temperature sensor	semiconductor (0-45°C)
Test available	FVC, VC, IVC, MVV, PRE-POST FVC, FEV1, FEV1/FVC%, DTPEF, FEV 0.5, FEV0.5/FVC%, FEV0.75, FEV0.75/FVC%, FEV2, FEV2/FVC%, FEV3, FEV3/FVC%, FEV6, FEV1/FEV6%, PEF, FEF25%, FEF50%, FEF75%, FEF25-75%, FEF75-85%, FET, Vext, ELA, EVOL, FIVC, FIV1, PIF, FIV1/FIVC%, FIF25%, FIF50%, FIF75%, R50, PIF, IRV, VC, IVC, IC, ERV, FEV1/VC%, TV, VE, RR, ti, te, ti/t-tot, TV/ti, MVV
Measured parameters	Up to 10000 tests

Memory capacity

Oximetry (on request)

Measurement method	Red and infrared absorption
SpO2 range	0-99%
SpO2 accuracy	± 2% between 70-99% SpO2
Average number of heart beats for the %SpO2 calculation	8 beats
Pulse Rate range	18-300 BPM
Pulse Rate accuracy	± 2BPM or 2% whichever is greater
Average interval for the calculation of cardiac pulse	8 seconds
Signal quality indication	0 - 8 segments on display
Test available	spot
Measured parameters	SpO2% min, max, average BPM min, max, average Test duration % Bradycardia Duration (<40 BPM) % Tachycardia Duration (>120 BPM) % of Time with SpO2 ≤ 90% (T90%, T89%)
Memory capacity	up to 300 hours oximetry

Certificates & Registrations

CE 0476	MED 9826
FDA 510 (k)	K 061712
Health Canada	71191 (class II), 75535 (class III)
CND code	Z12150102 (spiro) Z1203020408 (spiro + oxy)
GMDN code	46906 (spiro), 45607 (spiro + oxy)
Ministry of Health	1038066/R (spirometer) 1038086/R (spirometer + oximeter)

ITALY

MIR Head Office
Via del Maggiolino, 125
00155 Roma
Tel. +39 06 22 754 777
Fax +39 06 22 754 785
Mir.spirometry.com

USA

MIR USA, Inc.
5462 S. Westridge Drive
New Berlin, WI 53151
Phone +1 (262) 565-6797
Fax +1 (262) 364-2030

FRANCE

MIR Local Office
Jardin des Entreprises,
290, Chemin de Saint Dionisy
30980 LANGLADE (France)
Phone +33 (0)4 66 37 20 68
Fax +33 (0)4 84 25 14 32